

Turny Low Vehicle Swivel Seat

Turny Low vehicle electrical programmable swivel base that allows to rotate the seat and get it completely out of the vehicle. It allows the user to sit outside the vehicle and then be moved inside, comfortably and safely.



It is a **swivel base for low vehicles** that allows to take out the seat from the vehicle by 40 cm . Turny Low Vehicle is the most performant swivelling base currently marketed. It has unique features that makes it possible to use in car models where other swivel bases won't work



Turny Low provides a 4° inclination for improved comfort and pressure relief. The **Tilda function** optimizes the user's position during swivelling and extension.

Main features:

- Advanced technology swivel base
- Electric swivel and extension , easy to use
- Programmable stroke provides maximum leg and head space
- Adjustable comfort position
- 4° inclination for improved comfort
- 12 V POWER SUPPLY

Turny Low Vehicle: compatibility matrix

SEAT COMPATIBILITY				
BEV	BEV KIDS	COMPACT	RECARO	GS
✓	✓	✗	✓	✓ ¹

1. requires adapter plate

CARONY COMPATIBILITY			ACCESSORIES COMPATIBILITY		
CARONY	CARONY CLASSIC	CARONY GO	TILDA TILT	CARO-SLIDE	SLIDE RAILS
✗	✗	✗	✗	built in	✗

Turny Low Vehicle Swivel Seat



Address: BP19, 80490 Hallencourt, France



Tel : +33 (0)3 22 25 58 65 / Mobile : + 33 (0)6 45 72 46 90

contact@access-ability.eu
www.access-ability.eu

Turny Low Vehicle: models and dimensions

Turny Low Vehicle		Installation in vehicle	Length		Height	Swivel range	Movement, programming mode	Movement, comfort mode	Extension	Reverse extension (Re/Action [®])	Inclination	Weight capacity	Unit weight
			A	B									
103259	Turny Low Vehicle	right	61	47	11	0-113°	51	0-41	30	10	4°	140	65
103260	Turny Low Vehicle	left	61	47	11	0-113°	51	0-41	30	10	4°	140	65
103648	Application tool	right											
103649	Application tool	left											
103562	Universal mounting bracket*	front											
101865	Universal wireless remote												

All weights and dimensions are given in kilograms and centimetres, and for reference only. Weight capacity refers to the total weight of the seat, accessories and user.

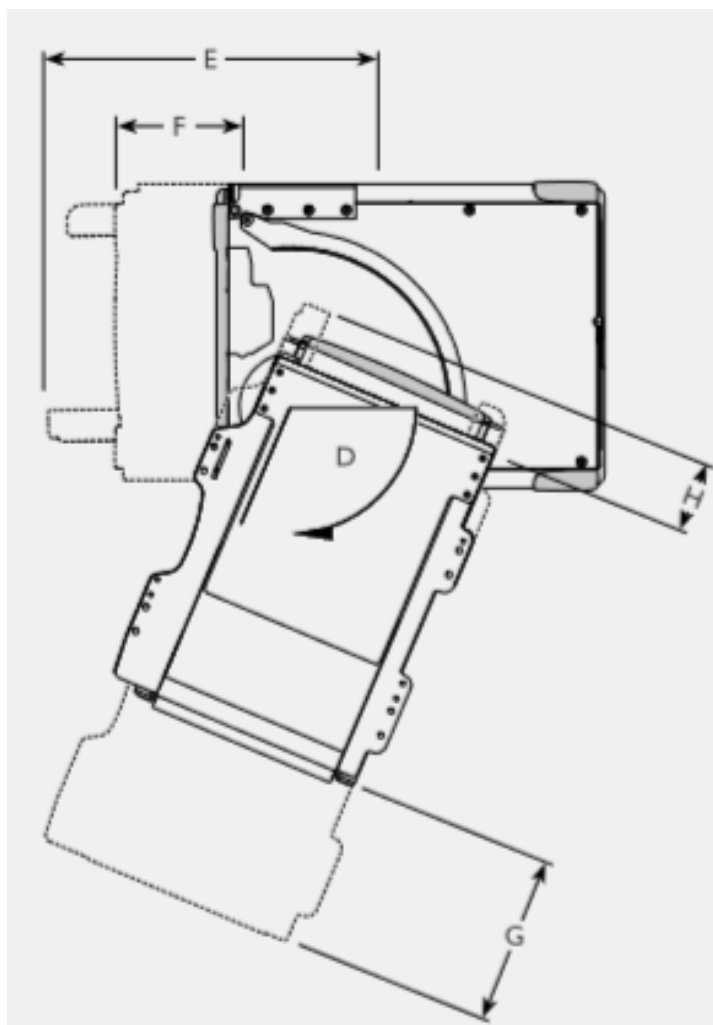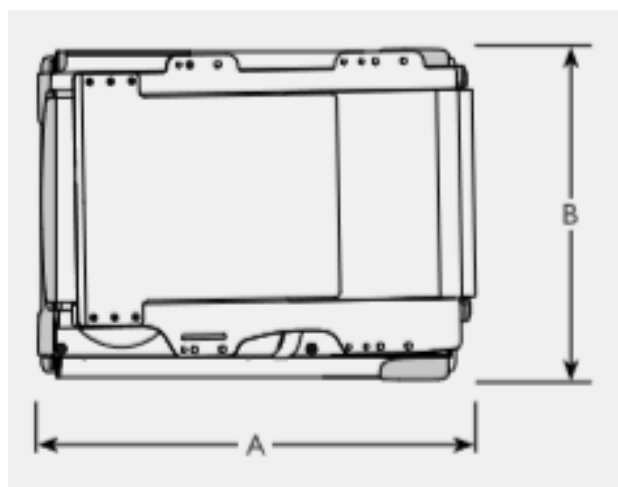

Turny Low Vehicle Swivel Seat

To guarantee your safety Autoadapt products are designed and tested according to current directives and standards.

Contact us for a free quotation!



Address : BP19, 80490 Hallencourt, France

Tel : +33 (0)3 22 25 58 65 / Mobile : + 33 (0)6 45 72 46 90

contact@access-ability.eu
www.access-ability.eu